



DIGITAL DISSOLUTION TESTER

DT LIGHT SERIES

www.erweka.com/products/dissolution-testers/usp-1-2-5-6-single-devices/dt-light/

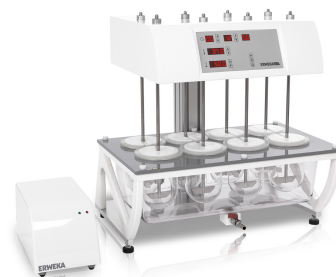
DT LIGHT SERIES

DIGITAL DISSOLUTION TESTER

The ERWEKA DT Light series offers the proven ERWEKA quality in a comprehensive package for a small budget, covering simple dissolution tests according to USP methods 1, 2, 5 and 6. The DTs are equipped with 6 or 8 test stations and a fixed drive head (high-head), which enables easy access to each 1000 ml vessel.

The shafts are easy to replace and the unique water bath made of molded PET is equipped with the tried and tested ERWEKA water outlet for easy cleaning. The external flow-through heater reduces the influence of vibrations and ensures a stable temperature. All accessories that come into contact with the dissolution sample are installed in Germany. Therefore, the DT Light series has the ERWEKA "Made in Germany" quality.

Together, these features ensure that the ERWEKA DT Light series is the perfect entry-level device into the world of dissolution testing.

**DT 126 LIGHT**

Perfect Entry-Level Dissolution Tester

- | High-head with 6 test stations
- | Dimensions (W x D x H): 510 mm x 450 mm x 660 mm
- | Weight: 30 kg

**DT 128 LIGHT**

Dissolution Testing With up to 8 Test Stations

- | High-head with 8 test stations
- | Dimensions (W x D x H): 642 mm x 482 mm x 680 mm
- | Weight: 38 kg

EASY START INTO THE WORLD OF DISSOLUTION TESTING

PRODUCT HIGHLIGHTS

**100% USP / EP / JP COMPLIANT**

Like all ERWEKA products, the DT Light series is 100% USP / EP / JP compliant

**HIGH-HEAD MODE**

The external flow-through heater reduces the influence of external vibrations and ensures a constant temperature

**USP METHOD 2**

Shafts with paddles for USP method 2 are included in the delivery

**EASY CLEANING**

Easy cleaning of the water bath and the set-up area

**EXTERNAL HEATER**

Like all ERWEKA products, the DT Light series is 100% USP / EP / JP compliant

**SIMPLE TO USE**

Simple control using a symbol keypad with LED display for temperature, rpm and run-time

DT LIGHT SERIES

TECHNICAL DATA

Weight	30 kg (DT 126 Light) 38 kg (DT128 Light)
Dimensions (H x W x D)	660 x 510 x 450 mm (DT 126 Light) 680 x 642 x 482 mm (DT128 Light)
Test stations	6 in 2 rows (DT 126 Light) 8 in 2 rows (DT 128 Light)
Sampling positions	High-head
USP methods	USP 1 / USP 2 / USP 5 / USP 6
Temperature	30 - 50 °C (± 0,2 °C)
Speed	20 - 250 rpm
Interfaces	1x RS-232, 2x USB

DT LIGHT SERIES

CONTACT US →

erweka.com/contact-services/contact-form/?prefilltext=contact_form.request&product=ER-D-S1-21

REQUEST A QUOTE →

erweka.com/contact-services/contact-form/?prefilltext=contact_form.quote_request&product=ER-D-S1-21

DT LIGHT SERIES PDF →

erweka.com/files/435519/dt-light-series.pdf

PRODUCT DATA SHEET →

erweka.com/files/462515/product-data-sheet-dt-light.pdf



ERWEKA
a **VERDER** company

www.erweka.com

www.erweka.com/products/dissolution-testers/usp-1-2-5-6-single-devices/dt-light/